
IDM

Instruments

Gelbo Flex Tester

Model No: G0002

Details: This apparatus is used to determine whether certain laminations of plastic withstand repetitive strain.

By attaching sample pieces of laminated plastic to the two circular clamping disks, via hose clamps, the specimen is twisted and turned to the specific movements. This process will show a visual result of pinholing and delamination.

Applications:

- Laminated plastic films.

Features:

- 45 rpm
- Stroke : 6"
- Twist : 400° Short stroke
- Twist : 440° Long stroke
- Counter for oscillations
- Totaliser

Options:

- Not applicable

Standards:


- ASTM F 392

Connections:

- Electrical: 220/240 VAC @ 50 HZ

Dimensions:

- **H:** 300mm • **W:** 1100mm • **D:** 350mm
- **Weight:** 45kg

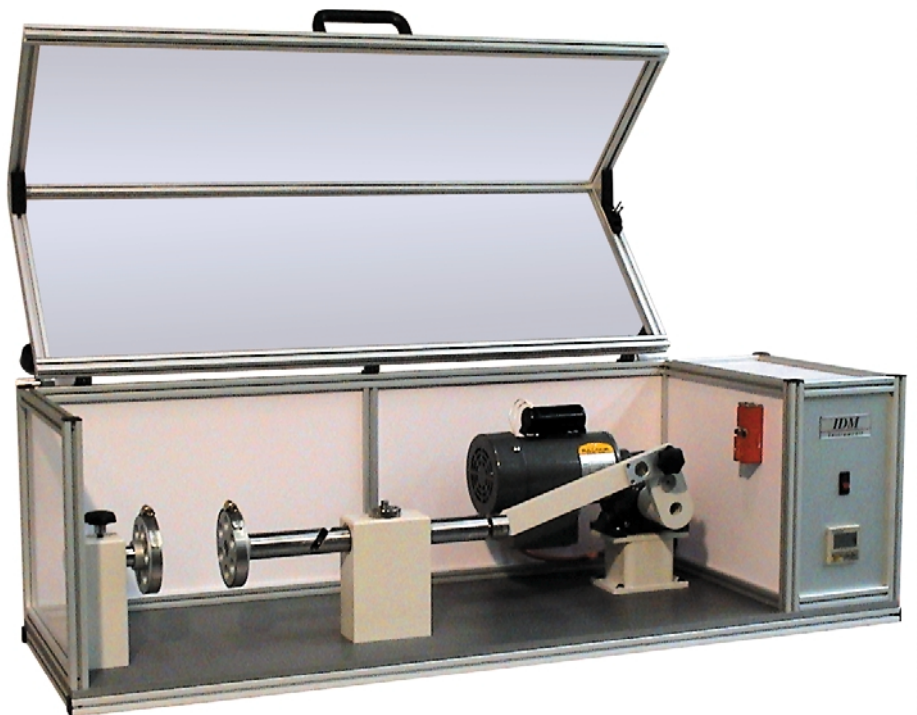

Representante Oficial para América Latina
REGMED INDÚSTRIA TÉCNICA DE PRECISÃO LTDA.

tel. +55 11 3601-5700 fax. +55 11 3601-5973 regmed@regmed.com.br www.regmed.com.br

All sizes and dimensions shown on this data sheet are averages only and images are as accurate as the medium allows.
Products and details shown may change without notice as IDM Instruments reserves the right to make alterations to its product range from time to time.